

## **TOOLS Prestained Protein Marker 10-245kDa**

**Cat.no.** TM-PM10245

**Storage:** Store at -20°C

**Product size:** 500 µl


### **Description**

TOOLS Pre-Stain Protein marker large is ready-to-use and pre-mixed with loading buffer for direct loading 5 µl on SDS-PAGE (Tris-glycine buffer).. This product does not need further dilution or addition of any reducing agent. After running the gel, the marker yields the following 12 discrete fragments (in KDa): 75(red), 25(green), 245, 180, 135, 100, 63, 48, 35, 20, 17, and 11(blue).

### **Recommendations**

1. Load 5µl/ gel well.
2. Do not boil the marker.

The product is for research only, not for diagnostic and clinical use.

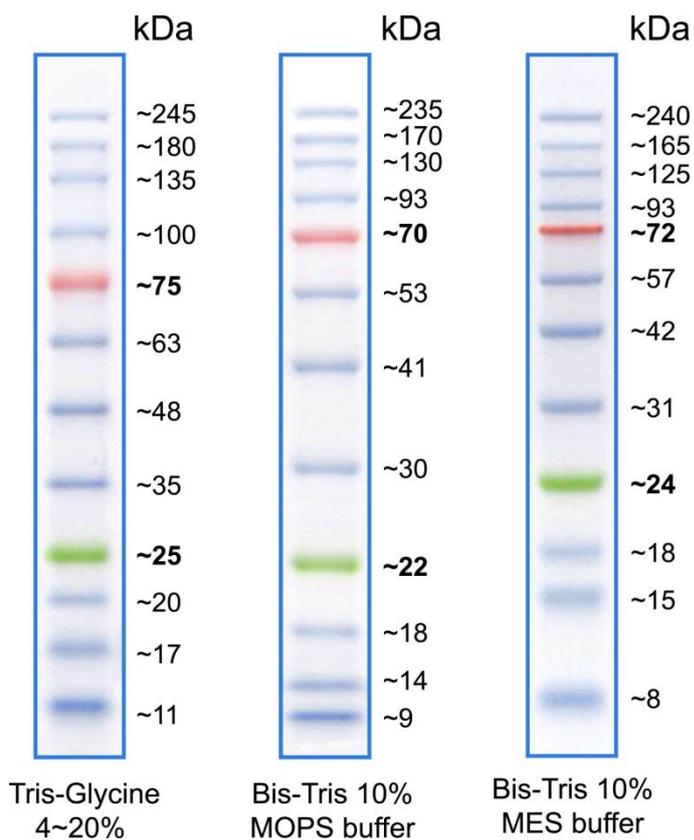

[www.tools-biotech.com](http://www.tools-biotech.com)

[info@tools-biotech.com](mailto:info@tools-biotech.com)

*Your best friend for research*

